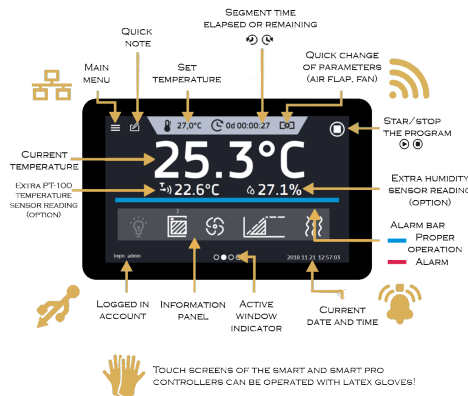




www.bunsen.es

DRYING OVEN



245L. DRYING OVEN WITH NATURAL AIR CONVECTION Ref. 8324 Plus

GENERAL CHARACTERISTICS

- * Interior and exterior construction in stainless steel AISI 304.
- * Solid door, fitted with a key. Excellent thermal insulation.
- * Protection against overtemperature class 3.1 according to DIN 12880.
 - * Large (7"), clear, full colour touch screen.
- * LAN, USB ports and WiFi for communication and data transfer.
 - * Multi-segment time and temperature programs.
 - * Loop function up to 255 times or endless.
 - * Adjustable start delay feature. Adjustable ramps.
- * Recording of min, average and max temperature value for each segment.
 - * Operating in temperature or time priority mode.
- * Temperature sensor fail alarm. Possibility of temperature calibration by the User.
- * Power failure control system (program continued after restoring power).
 - * Overview of data in tabular and graphic form.
 - * Visual and sound alarm.
 - * Administrating functions for management.
- * Password protected log-in. Internal memory for programs and data storage.
- * Operating with gloves on. Event registry with user notifications.




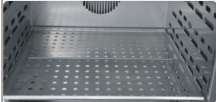

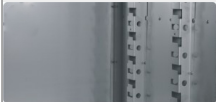




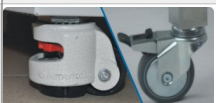
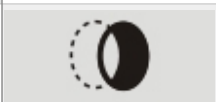
BUNSEN SA. 30 Rio Alberche St - 28970 Humanes de Madrid, Madrid, SPAIN
Tlf. +34 916 113 584 - Fax. +34 916 128 254
info@bunsen.es - www.bunsen.es



- * Alarm Bar – quick visual information about device status.
- * Quick Note – user can save messages (50 characters) in Smart PRO controller memory.
 - * Quick Change of parameters: temperature, humidity, time, air flap and fan.
- * Trays made of stainless steel wire; 3 tray included, maximum capacity; 10. Max shelf workload 25 kg. Consult accessories.

Manufactured according to CE directives.

| Ref. | Capacity | Dimensions W x H x D | Useful measures W x H x D | Weight | Temperature | Nominal power |
|-------|----------|----------------------|---------------------------|---------|-------------|---------------|
| 8324P | 245L. | 810 x 1200 x 770mm. | 600 x 800 x 510mm. | 126 Kg. | Up to 300°C | 3100 |

| Ref | Accessories | |
|---------|---|--|
| 8300-02 |  | Stainless steel wire shelf INOX |
| 8300-03 |  | Perforated shelf |
| 8300-04 |  | Reinforced shelf |
| 8300-05 |  | Reinforced version |
| 8300-06 |  | Stainless steel cuvettes |
| 8300-07 |  | Table with wheels |
| 8300-08 |  | Additional Pt 100 temperature sensor |
| 8300-09 |  | HEPA-fresh air filter |
| 8300-10 |  | Wheels |
| 8300-11 |  | Non-standard access port for external sensor |



BUNSEN SA. 30 Rio Alberche St - 28970 Humanes de Madrid, Madrid, SPAIN
 Tlf. +34 916 113 584 - Fax. +34 916 128 254
info@bunsen.es - www.bunsen.es

