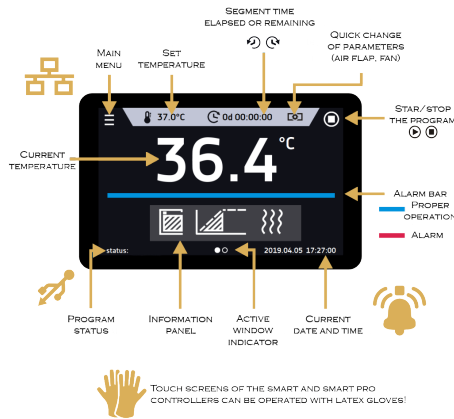




[www.bunsen.es](http://www.bunsen.es)



## DRYING OVEN

### 53L. DRYING OVEN WITH NATURAL AIR CONVECTION Ref. 8302

#### GENERAL CHARACTERISTICS

- \* Interior and exterior construction in stainless steel AISI 304.
- \* Solid door, fitted with a key. Excellent thermal insulation.
- \* Protection against overtemperature class 2.0 according to DIN 12880.
  - \* 4,3", clear, full colour touch screen.
  - \* USB and LAN ports for data transfer
- \* Multi-segment time and temperature programs. Loop function up to 255 times or endless.
  - \* Adjustable start delay feature, adjustable ramps.
- \* Recording of min, average and max temperature value for each segment overview of set and current parameters while operating.
  - \* Operating in temperature or time priority mode.
- \* Fan speed control in the range of 0 ... 100% or 10% ... 100% (depending on the chamber capacity).
  - \* Temperature sensor fail alarm.
- \* Possibility of temperature calibration by the User.
- \* Power failure control system (program continued after restoring power).
  - \* Visual and sound alarm.



BUNSEN SA. 30 Rio Alberche St - 28970 Humanes de Madrid, Madrid, SPAIN  
Tlf. +34 916 113 584 - Fax. +34 916 128 254  
[info@bunsen.es](mailto:info@bunsen.es) - [www.bunsen.es](http://www.bunsen.es)




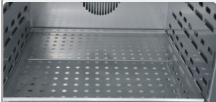

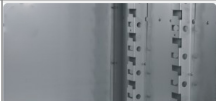




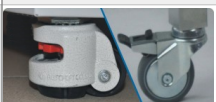
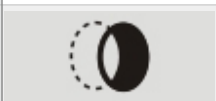
\* Quick change of program parameters.

\* Alarm Bar.

\* Shelves standard included 2, possibility of including 5. Max shelf workload 25Kg.

Manufactured according to CE directives.

Ref.	Capacity	Dimensions W x H x D	Useful measures (W x H x D)	Weight	Temperature	Nominal power
8302	53L.	590 x 700 x 620mm.	400 x 390 x 360mm.	50 Kg.	Up to 300°C	1700

Ref	Accessories	
8300-02		Stainless steel wire shelf INOX
8300-03		Perforated shelf
8300-04		Reinforced shelf
8300-05		Reinforced version
8300-06	1, 2, 3, ... 	Door openings counter
8300-07		Table with wheels
8300-08		Additional Pt 100 temperature sensor
8300-09		HEPA-fresh air filter
8300-10		Wheels
8300-11		Non-standard access port for external sensor



BUNSEN SA. 30 Rio Alberche St - 28970 Humanes de Madrid, Madrid, SPAIN

Tlf. +34 916 113 584 - Fax. +34 916 128 254

[info@bunsen.es](mailto:info@bunsen.es) - [www.bunsen.es](http://www.bunsen.es)

